Hoffen

BALANÇA DE ANÁLISE CORPORAL BODY COMPOSITION ANALYSER SCALE

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English version: page 19-35

Manual version: 11.2024

NOTE: Item colours may vary slightly from those shown.

BEFORE USING THE APPLIANCE, PLEASE READ THE INSTRUCTION MANUAL CAREFULLY, ESPECIALLY THE SAFETY INSTRUCTIONS, AND KEEP IT FOR LATER USE. IF YOU GIVE THE APPLIANCE TO A THIRD PARTY, PLEASE ALSO GIVE THEM THIS INSTRUCTION MANUAL.

1. INTRODUCTION

Thank you for choosing us. This appliance will certainly make your day-to-day life easier. Please read this instruction manual carefully before putting it into operation so that you can get the most out of it, and keep it for future reference. Keep this manual for reference when using the product.

2. TECHNICAL INFORMATION

Material: tempered glass, plastic and stainless steel

Dimensions: 26x28 cm

Batteries: 2 x 1.5V AAA

Capacity: max 180 kg

Units: kg / LB

Memory: up to 8 users

Analysis: Weight, BMI, body fat, muscle, water, visceral fat, bone mass, metabolism

Indicator: Overload indicator, Low battery indicator

APP: OKOK

3. GENERAL INFORMATION

This scale is a high-precision instrument for measuring and analysing body composition data.

You can transmit your body composition data to your mobile phone using Bluetooth technology. You can graph your weight and other data, from which you can follow the trend at any time. Download the app to your mobile phone and start monitoring your weight and body shape. The product image is purely indicative. Appearance specifications must be respected.

4. SAFETY WARNINGS

- To guarantee the accuracy of the data, it is

recommended that each measurement be taken at the same time of day.

- Do not use in conditions of strong electromagnetic field interference.
- If the scale will not be used for a long time, when you make the first weighing, touch the scale so that the display shows zero and let the scale switch off automatically. You can then touch it again to start weighing.
- To ensure the accuracy of the measurement, take off your shoes, socks and other heavy objects, keep your hands and feet clean, then step onto the scale with bare feet and hold the handle with both hands.
- Only use the appliance on stable, flat surfaces that are protected from impacts, and never use it on carpets or rugs. Placing the scale on a hard, flat surface will ensure greater precision and consistency.
- Do not use the appliance with damp or wet feet.
- The platform is slippery when wet. Keep it dry for your safety.

- Climb onto the weighing platform with care. Do not climb with only one foot onto the edge of the scale at the risk of tipping it over; always lean on the centre of the scale and adjust your feet to touch the 4 electrodes once you are stable. Stay steady and balanced to avoid slipping.
- It is not recommended to install the scale in the following places: inclined or unstable surfaces; places that suffer strong interference from electromagnetic fields; places with vibration, rotation and mobile equipment; places with draughts and large temperature variations.
- Do not use the scale in humid or overheated environments (temperature range 5° to 35°C).
- This device is a highly sensitive precision electronic instrument. Keep it away from high-frequency devices as they may cause interference. If you detect malfunctions due to interference, move the device away from that area or switch off the device that is affecting its operation.
- Measured results may be misleading after intensive physical exercise, excessive dieting or

under extreme conditions of dehydration.

- With regard to analysing body composition, the measurement data may not be accurate for the users listed below:
 - Children under the age of 18 or adults over the age of 70;
 - Professional sportspeople, such as bodybuilders.
- The use of this product is not recommended for the following people:
 - Pregnant women;
 - People with a pacemaker, electrocardiograph or any other device implanted in the body or life support device;
 - People with symptoms of oedema or who are undergoing dialysis are prohibited from using this device.
- To avoid damaging the device, make sure that no liquid comes into contact with the electrical parts. Do not use the appliance outdoors and keep it away from direct or indirect sources of water. If any liquid spills on the electrical parts, switch off the appliance and contact the shop where you

purchased it or a qualified technician for repair.

- The use of this appliance by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience and knowledge in handling it is not recommended, unless they are given supervision and instructions on how to use the appliance by the person responsible for their safety when using it.
- Children are unaware of the damage that can be caused by electrical appliances and must always be supervised to ensure that they do not play with the appliance.
- This appliance is not a toy and must therefore be kept out of the reach of children.
- After unpacking the appliance, make sure that it is not damaged and that no parts are missing. If in doubt, contact the shop where you purchased the appliance. All contents of the packaging, such as plastic bags, styrofoam, tape, etc., must be kept out of the reach of children as they can be dangerous and cause suffocation.
- Never dismantle or attempt to repair the appliance

in any way. Failure to comply with this safety standard could jeopardise your own safety and invalidate the terms of the guarantee. Only suitably qualified technical personnel may carry out any type of repair on this appliance.

- The appliance has no user-serviceable components. Any type of repair must be carried out by suitably qualified technical personnel. Ensure that only new and recommended parts are used.
- Always keep the scale in a horizontal position.
- Never clean the device with abrasive cleaning detergents or scouring pads, alcohol or solvents.
- The scale is an electronic instrument and should never be submerged in water.


BATTERIES

- **WARNING:** this appliance requires 2 x AAA batteries (2x1.5V). Batteries can cause serious burns if swallowed, potentially resulting in internal bleeding or death. Keep out of sight and reach of small children and ensure the battery compartment is always properly closed.

- When inserting batteries into the appliance, make sure your hands are dry.
- Insert 2 AAA batteries into the battery compartment, respecting the polarity.
- Batteries may only be replaced with batteries of the same type or equivalent.
- Remove the batteries from the appliance when not in use for a long period of time and store them in a cool, dry place.
- DANGER OF EXPLOSION! Batteries must not be recharged or reactivated by other means, dismantled, thrown into a fire or short-circuited.
- Remove exhausted batteries from the appliance immediately, as they may leak and cause damage.
- WARNING: Batteries can be life-threatening if swallowed. Therefore, keep batteries out of the reach of children. If a battery has been swallowed, medical attention should be sought immediately.
- Do not expose batteries to extreme conditions, do not solder them, do not place them on heat sources or expose them to direct sunlight, as these conditions will damage them.

- If the batteries leak acid, avoid contact with the skin, eyes and mucous membranes. In the event of contact with acid, immediately rinse the affected area with plenty of water and seek medical attention.
- Check the batteries regularly. Leaking batteries can damage the appliance.
- Wear protective gloves if the batteries leak liquid. Clean the battery compartment and contacts with a dry cloth.
- After use, the batteries must be disposed of properly.

5. APPLIANCE DESCRIPTION

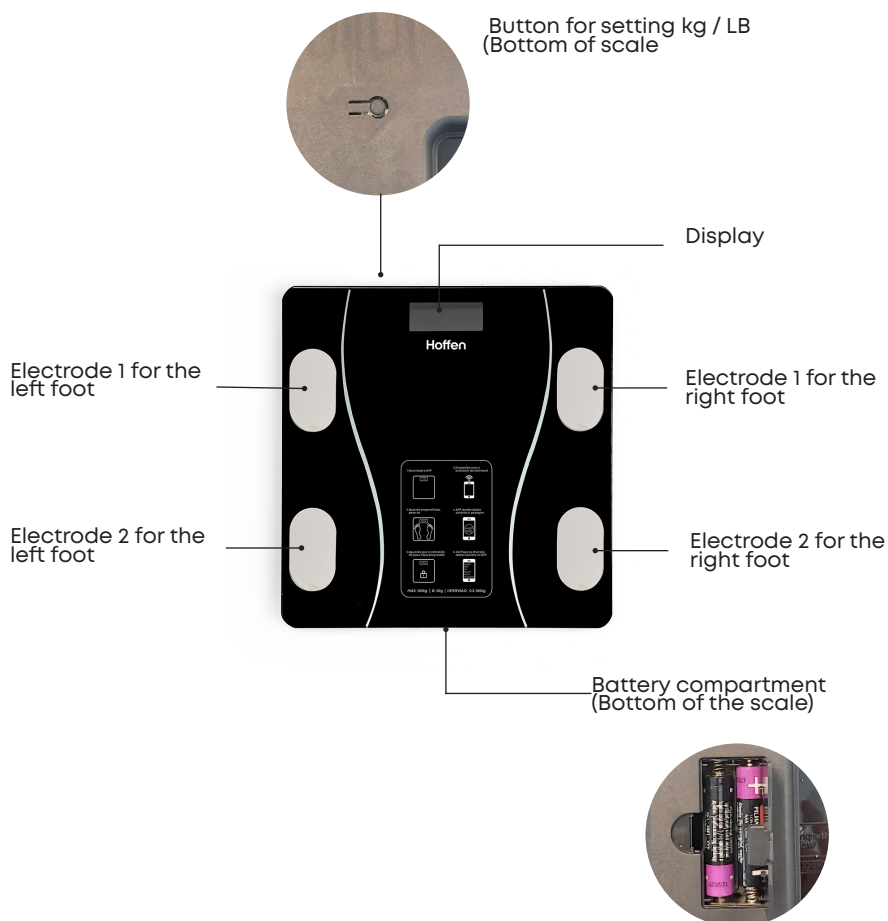6. INSTALLATION AND BATTERY REPLACEMENT

- The scale requires 2 x 1.5V AAA batteries, which must be installed as follows:

First use: Remove the insulation film before use (Note: if batteries are included). If batteries are not included, open the battery compartment cover, install 2 batteries respecting polarity and close the compartment

properly.

Replacing the batteries: Open the battery compartment cover, remove the old batteries and then install the new ones, respecting the polarity.

Note: After installing the batteries, close the compartment lid. Please use the batteries specified for this appliance. Remove the batteries if the appliance will not be used for a long time (more than three months).

7. DOWNLOAD, APP INSTALLATION AND REGISTRATION

OKOK INTERNATIONAL APP

This application is supported by Android 4.3 and above and selected Android devices with Bluetooth 4.0 and above. The application can store data from up to 8 people, so you can use this scale with your family and friends.

- Search for 'OKOK' in the App Store or Google Play to download the application or, more conveniently, scan the QR code below.

Download now:
OKOK International

ANDROID APP ON
Google play

Available on the
App Store

OKOK Healthcare

Tips: The application (APP) is updated continuously, so it may differ from what is shown in the manual. For final and accurate information, check on your mobile phone.

INSTALLATION AND CONNECTION

- The device must be connected to the application before it is used for the first time.
- The scale automatically connects to your mobile phone and the app.
- Make sure that Bluetooth is activated on your mobile phone before you start using the scale.
 1. Switch on Bluetooth and use your device to scan the QR code or search for OKOK International in the App Store/Google Play to install the app.
 2. The mobile phone or tablet connects to the network (WIFI, 3G, 4G, 5G).
 3. Pair the device, register and complete the user information (enter name, gender, height, date of birth and preferred weight). Click 'Next' to confirm your personal information.

4. Weigh and synchronise. Get on the scales and weigh yourself. Wait for the device to calculate your measurements.
 5. Visualise the result.
- The device displays body weight, BMI, fat, muscle, water, visceral fat, bone mass and metabolism. The device automatically saves your data.
 - The device displays your current weight, generates reports and gives health advice based on your measurements. The weight trend includes a weekly, monthly and quarterly report, indicating weight changes through curves. Account management includes registration, logging in, changing the password, forgetting the password, connecting to social networks and much more. User management includes information about the user, such as nickname, birth date, height, gender, etc.

8. CLEANING

- Do not use abrasive cleaning detergents, scouring pads or metal utensils to clean the appliance.
- Wipe the outside of the scale with a dry or slightly damp cloth and then dry it with a clean, dry cloth.
- Never immerse the appliance in water or any other liquid.

9. REPAIRS

- This appliance has no user-serviceable parts. Do not attempt to repair the appliance yourself.
- If the power cable is damaged, it must be replaced by the manufacturer, an authorised service agent or a qualified professional.

10. DISPOSAL OF WASTE EQUIPMENT

Waste electrical devices are made of recyclable materials. Do not dispose of them with household waste as they may contain substances harmful to human health and the environment! Please actively support cost-effective management of natural resources and protect the environment by taking the waste appliance to a collection point for recyclable waste electrical and electronic equipment.

11. DECLARATION OF CONFORMITY EU

This appliance complies with the requirements of the Electromagnetic Compatibility Directive (EMC) and the Low Voltage Directive (LVD), therefore it was marked with the CE label and provided with a declaration of conformity with European standards.

12. WARRANTY

The product warranty is valid for three years. In the event of malfunction, the product must be returned to the store where it was purchased, along with the proof of purchase. The product should be complete and placed in its original packaging.

13. SYMBOLS

This symbol indicates that the user manual is included in the product packaging and must be read before use.

	The „Green Dot” trademark located on the packaging means that the entity, for whom the product was manufactured, made financial contributions in construction and operation of the national recovery and recycling system for packaging waste in accordance with the principles resulting from the laws of Portugal and the European Union concerning packaging and packaging waste.
	This symbol means that the product must not be disposed of with other household waste. Portuguese law prohibits, under penalty of a fine, combining waste electrical and electronic equipment with other wastes. It is possible to prevent potential negative consequences for the environment and human health, which could otherwise arise from inappropriate handling of waste electronic equipment, if the product is disposed of properly.
	The ‘CE’ symbol is the guarantee of compliance with harmonised European standards, on a voluntary basis, reflecting the essential requirements in the technical specifications and is intended to ensure that the appliance complies with the requirements that provide a high level of protection for the health and safety of people, domestic animals and property, while guaranteeing the functioning of the internal market.
	This product meets the requirements of Directive 2011/65/EU, referred to as RoHS directive. The purpose of the Directive is the approximation of the laws of the member states concerning the restriction of use of hazardous substances in electrical and electronic equipment, and the contribution towards the protection of human health and the environmentally sound recovery and disposal of waste electrical and electronic equipment. RoHS Directive applies from 3 January 2013. New electrical and electronic equipment put on the market shall not contain lead, mercury, cadmium, hexavalent chromium, polybrominated biphenyls or polybrominated diphenyl ethers.

	Batteries should not be disposed of with household waste. The disposal of batteries is regulated by national and European laws. Batteries must be collected separately and recycled because they contain heavy metals (e.g. lead, mercury and cadmium) that can cause serious environmental and health problems. As a consumer, you have a legal duty to separate your used batteries from your household rubbish and take them to a collection point for recycling.
	This selective disposal symbol indicates that packaging components made of cardboard or paper should be placed in the blue-coloured recycling bin for these materials.

14. EQUIPMENT PRODUCED FOR:

Joinco, Import e Export, Lda.
Ed. ARCIS, Rua Ivone Silva 6 15º piso
1050-124 Lisbon
PORTUGAL
VAT 507 191 765
www.joinco.pt
Made in China

15. EUROPEAN DECLARATION OF CONFORMITY

Joinco, Import & Export, Ltd.,
Ed. Arcis, Rua Ivone Silva nº6, 15º piso
1050-124 Lisboa, Portugal
declares that the equipment:

Description: Body composition analyser scale

Model: iBDS-H209

Brand: Hoffen

Lot no.: PD/559/2024

APP: OKOK

is in compliance with Directive 2014/53/EU of the European Parliament and of the Council of 16 April 2014.

The full text of the declaration of conformity is available at the following internet address: **www.joinco.pt. joinco.pt**

There are no known restrictions on putting into service or requirements for authorisation of use in any EU country.

Operating frequency band: 2402-2480 MHz

Maximum transmitted power (EIRP): 12 dbm

MANUFACTURED BY:

Joinco, Import e Export, Ltd.

Ed. ARCIS, Rua Ivone Silva 6 15º piso

1050-124 Lisbon

PORTUGAL

VAT 507 191 765

www.joinco.pt

Hoffen
ALL FOR BETTER LIFE